



TEL:4006-871-227 Web:www.ybio.net Email:shybio@126.com

YBB790Hu01 50μg

Recombinant Sialic Acid Binding Ig Like Lectin 8 (SIGLEC8)

Organism Species: Homo sapiens (Human)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

10th Edition (Revised in Jan, 2014)

[PROPERTIES]

Residues: Gly70~Ser336

Tags: Two N-terminal Tags, His-tag and GST-tag

Accession: Q9NYZ4

Host: *E. coli*

Subcellular Location: Membrane; Single-pass type

I membrane protein.

Purity: >95%

Endotoxin Level: <1.0EU per 1 μg

(determined by the LAL method).

Formulation: Supplied as lyophilized form in 20mM Tris, 150mM NaCl, pH8.0, containing 1mM EDTA, 1mM DTT, 0.01% sarcosyl, 5% trehalose, and preservative.

Predicted isoelectric point: 9.2


Predicted Molecular Mass: 59.0kDa

Applications: SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

[USAGE]

Reconstitute in sterile ddH₂O.

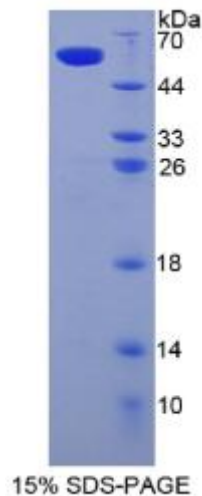



TEL:4006-871-227 Web:www.ybio.net Email:shybio@126.com

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCES]

The sequence of the target protein is listed below.

G DRPYQDAPVA TNNPDREVQA ETQGRFQLLG DIW SNDCSLS IRDARKRDKG
SYFFRLERGS MKWSYKSQLN YKTKQLSVFV TALTHRPDIL ILGTLESGHS RNLTCSVPWA
CKQGTPPMIS WIGASVSSPG PTTARSSVLT LTPKPQDHGT SLTCQVTLPG TGVTTTSTVR
LDVSYPWNLT TMTVFQGDAT ASTALGNGSS LSVLEGQSLR LVCAVNSNPP ARLSWTRGSL
TLCPSRSSNP GLLLEPRVHV RDEGEFTCRA QNAQGS

[REFERENCES]

1. Varela I., *et al.* (2011) Nature 469:539-542.
2. Na H.J., *et al.* (2012) J. Immunol. Methods 383:39-46.
3. Na H.J., *et al.* (2012) Cytokine 57:169-174.
4. Hudson S.A., *et al.* (2011) J. Clin. Immunol. 31:1045-1053.