



TEL:4006-871-227 Web:www.ybio.net Email:shybio@126.com

YBJ103Mu01 10 μ g

Recombinant Filaggrin (FLG)

Organism Species: *Mus musculus* (Mouse)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

[PROPERTIES]

Residues: Asp119~Tyr336 (Accession # P11088), with two N-terminal Tags, His-tag and GST-tag.

Host: *E. coli*

Purity: >95%

Endotoxin Level: <1.0EU per 1 μ g (determined by the LAL method).

Formulation: Supplied as lyophilized form in PBS, pH7.4, containing 5% sucrose, 0.01% sarcosyl.

Predicted isoelectric point: 6.7 Predicted Molecular


Mass: 55.8kDa

Applications: SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

[USAGE]

Reconstitute in sterile PBS, pH7.2-pH7.4.

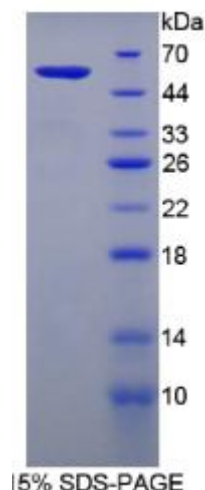

5% SDS-PAGE



[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCES]

The target protein is fused with two N-terminal Tags, His-tag and GST-tag, its sequence is listed below.

MSPILGYWKI KGLVQPTRLL LEYLEEKYEE HLYERDEGDK WRNKKFELGL EFPNLPYYID
GDVKLTQSMA IIRYIADKHN MLGGCPKERA EISMLEGAVL DIRYGVSRIA YSKDFETLKV
DFLSKLPEML KMFEDRLCHK TYLNGDHVTH PDFMLYDALD VVLYMDPMCL DAFPKLVCFK
KRIEAIQID KYLKSSKYIA WPLQGQWATF GGGDHPPKSD GSTSGSGHHH HHHSAGLVPR
GSTA IGMKET AAA KFERQH M DSPDLGTLEV LFQ GPLGS EF-DF SEGQ AVG AHR
QSGAG QRH EQ RSSRGQ HGSG YYY EQEHSEE ESDS QHQH GH QHEQQ RGHQ H
QHQHQHEHEQ PESGHRQ QQS SG RGHQGAHQ EQ GRDSARPR GSNQG HSSSR
HQ A D S P R V S A R S G S G G R G Q S P D A S G R S S N R R D R P R Q P S P S Q S S D S Q V H S G
VQVEAQRGQS SSANRRAGSS SGSGVQGAAA SGQGGYESIF TAKHLDFNQS HSYYYY

[REFERENCES]

1. Rothnagel J.A., *et al.* (1987) J. Biol. Chem. 262:15643-15648.
2. Eddleston J., *et al.* (1999) Genomics 56:149-159.



TEL:4006-871-227 Web:www.ybio.net Email:shybio@126.com

3. Yang A., *et al.* (1999) Nature 398:714-718.

4. Presland R.B., *et al.* (2000) J. Invest. Dermatol. 115:1072-1081.